Sonifer

TERMÓMETRO DIGITAL – MANUAL DE INSTRUCCIONES
DIGITAL THERMOMETER - INSTRUCTION MANUAL
THERMOMETRE DIGITAL - MANUEL D'INSTRUCTIONS
TERMÔMETRO DIGITAL - MANUAL DE INSTRUÇÕES

TEM 6000

Sonifer, S.A.
Avenida de Santiago, 86
30007 Murcia
España
E-mail: sonifer@sonifer.es

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ENGLISH

WHAT IS A "NORMAL" TEMPERATURE?

Although the generally accepted "normal" reading is 37.0°C temperatura readings can vary from 36.1 °C/97.0 to 37.2°C and still may be considered as "normal" . Variations in temperature can be attributed to activities such as exercise , smoking , eating and drinking . Body temperature is lower in the morning than in the afternoon . Other variations are due to the body where temperature is measured . While oral temperature follows the guidelines previously described , rectal temperature is generally 0.5°C of higher, and axillary(under the arm)temperature will be 0.5°C of lower.

SPECIFICATIONS

Measurement range: 32 °C ---42. 9°C

Low temperature display: Temperature<32 °C display LOC

High tern perature display: Temperature > 43 °C display H°C

Accuracy: 35---42 °C +/- 0.1 ,(within the ambient temperature range from 18°C to 28°C).

Display resolution: 0.1°C

Normal beep sound: 10 long beep sounds in 10 seconds

Fever beep sound(>3 7.8°C: 30 short beep sounds in 10 seconds

Auto power off: 8minute +/- 2minute

Battery: AG3 or LR41 ; Life of battery: 2years or 2500times

Dimension: 125mm X 19mm X 11mm

Unit Weight: approx X 9.8g (including battery)


Operating condition: Temperature 1 0°C---40 °C

Humidity 30%---85%RH

Storage condition: Temperature -10°C ---60°C

Humidity 25%---90%RH

HOW TO TAKE YOUR TEMPERATURE

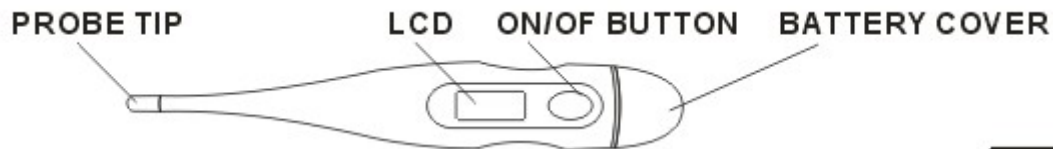1. Press ON/OFF button to activate . The unit would beep and display
This LCD display test will go on for around 2 seconds.
2. Then, the temperature taken last time would be displayed with a "M" mark in the right corner.
3. When L and flashing °C displays , the thermometer is ready for temperature measuring . But if the room temperature is lower than 32 .0°C the room temperature would be displayed instead of "L".
4. The beep will sound when temperature measuring is complete .
5. Press ON/OFF button of the unit. (The unit will automatically power off after approximately eight minutes after use.)

TAKING AN ORAL TEMPERATURE: Place the probe tip under the tongue as near as possible to the heat pocket. Approximate measuring time: 1 minute .

TAKING AN AXILLARY TEMPERATURE (UNDER THE ARM): Wipe the underarm with a dry towel and place the probe tip under the arm so the tip is touching the skin the thermometer perpendicular to the body. Position the arm across the chest so the probe tip is well covered by the arm. Approximate measuring time :2-3minutes.

Note: axillary temperatures are generally 0.5°C/1 .0°F lower than oral readings.

TAKING A RECTAL TEMPERATURE: Apply a water-soluble lubricant to the probe cover. Gently insert the probe (MAXIMUM 1/2") into the rectum . Approximate measuring time :1 minute.

Note: Rectal temperature is generally 0.5°C/1 .0 of higher than oral readings.

REPLACING THE BATTERY

The battery in your digital thermometer needs to be replaced when "🔋" appears in the lower right corner of the display. Replace the battery with a 1.5V, AG3 alkaline or LR41 silver oxide type or equivalent.

To replace the battery:


1. Pull the battery cover away from the unit.

2. Pull up the battery set about 18mm.

3. Use suitable tool to remove the old battery and change a new one.

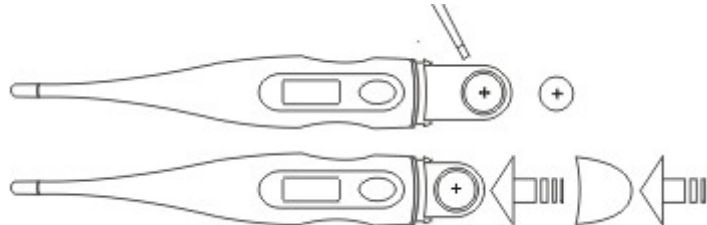CLEANING INSTRUCTIONS

Before and after each use, clean the thermometer with a soft cloth and isopropyl alcohol diluted with water, or wash with COLD soapy water. The unit can be submerged in water or alcohol (for waterproof model only). Do not soak for an extended period of time or sterilize by boiling method, gas or steam autoclaves.

GENERAL CARE

1. The thermometer should only be used under the supervision of an adult.
2. Do not walk, run or talk during temperature taking
3. Store the unit in the protective case when not in use.
4. Do not store the unit where it will be exposed to direct sunlight, dust or humidity. Avoid extreme temperatures.
5. Dropping or subjecting your thermometer to strong shocks should be Avoided
6. Do not attempt to disassemble the unit, except for replacing the Battery.

Device classification: Class BF

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

Kangfu Medical Equipment Factory
No. 380 Ningkang East Road, Yueqing, Zhejiang 325600, China

ROVAL MED SRL
Industiilor, nr. 2M, Bucharest, 032895-Romania

2010

Lot: KFT-04201264-66

CE 0197

Instruction manual: 30/07/2012

Disposal of old electrical appliances.

The European directive 2002/96/EC on Waste Electrical and Electronic Equipment (WEEE), requires that old household electrical appliances must not be disposed of in the normal unsorted municipal waste stream. Old appliances must be collected separately in order to optimize the recovery and recycling of the materials they contain, and reduce the impact on human health and the environment.

The crossed out “wheeled bin” symbol on the product reminds you of your obligation, that when you dispose of the appliance, it must be separately collected.

Consumers should contact their local authority or retailer for information concerning the correct disposal of their old appliance.

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